

3D MaxiPol Set



Ref.No 313240

Polieren von Prothesen- und 3D-Kunststoffen
am Poliermotor

Polishing of acrylic dentures and 3D-resins
on the polishing lathe

Polissage des prothèses acryliques et résines 3D
sur le moteur de polissage



Politur am Poliermotor

on the polishing lathe

sur le moteur de polissage

1

Vorpolitur

pre-polishing

pré polissage




Einlagebürste
4280

+



Polistar
Lintygrau 90g

→

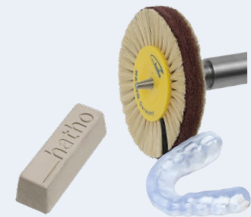

2

Glanzpolitur

shine polishing

poli brillant



Mikrofaserleder
520 100/7

+



Polistar
Beige 85g

→



3

Hochglanzpolitur

high shine polishing

lustrage



Baumwollfaden
152 80 BW

+



Polistar
Emulsion 50ml

→



Info

4280 = Ø 80mm = 2.800 rpm
520 100/7 = Ø 100mm = 1.400 rpm
152 80 BW = Ø 80mm = 2.800 rpm

3D MaxiPol Set

Inhalt/contents/contenu

2x 4280, 1x Polistar Lintygrau 90g, 1x 520 100/7, 1x Polistar Beige 85g,
1x 152 80 BW, 1x Polistar Emulsion 50ml